



Field Installation and Usage Instructions

Part No. TMA29-15

Paper Roll Holder Assembly

For All Model TMM4 Chairs



5516 S.W. 1st Lane – Ocala, FL 34474-9307
Phone: 1-800-237-3377 - Fax: 1-352-854-9544
www.transmotionmedical.com

Document No.: IM TMA29-15
Revision: C

TABLE OF CONTENTS

TABLE OF CONTENTS.....	2
IMPORTANT NOTES ABOUT THIS MANUAL.....	3
SPECIAL NOTES - SIGNAL WORDS.....	3
TOOLS REQUIRED.....	4
PARTS INCLUDED IN TMA29-15 KIT.....	4
PREPARING THE BACKREST FOR THE PAPER ROLL HOLDER.....	4
INSTALLING THE PAPER ROLL HOLDER ASSEMBLY.....	6
USING THE PAPER ROLL HOLDER ASSEMBLY.....	8
CLEANING INSTRUCTIONS.....	12

IMPORTANT NOTES ABOUT THIS MANUAL

This Field Installation and Usage Instruction document is designed to assist you with installation of the TransMotion Medical Product. Carefully read this document thoroughly before using the equipment or doing service/maintenance on it. If you are unable to understand the WARNINGS and CAUTIONS and instructions, contact TransMotion Medical customer service before attempting to operate or service the equipment otherwise injury or damage may result.

To assure safe operation of this device, it is essential that methods and procedures be established for educating and training staff on the safe and effective operation of the product.

TransMotion Medical reserves the right to change specifications without notice.

The information contained in this document is subject to change without notice.

Specifications listed are nominal and may vary slightly from unit to unit due to tolerance variations or power supply variations.

Photographs and illustrations contained within this document may not depict exactly the model you have. This document can cover many closely related models.

SPECIAL NOTES - SIGNAL WORDS

Signal words are used in this manual and apply to hazards which could result in injury or property damage. The following is a definition of those signal words as used in this document:

DANGER

Indicates an imminently hazardous situation which, if not avoided, will result in death or serious injury.

WARNING

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices or potential property damage hazards.

NOTICE

Provides important information, makes special instructions clearer, or provides service personnel information to make maintenance easier.

TOOLS REQUIRED

The tools needed to install the Paper Roll Holder are as follows:

7/16" Wrench
Pencil
5/32" Allen Wrench
1/4" Drill Bit
Power Drill




PARTS INCLUDED IN TMA29-15 KIT

Verify that you have all the parts shown. If not, contact TMM Service at 1-800-237-3377.


- (1) Paper Roll Bracket Assembly
- (1) TMM-559-11 – Paper Roll Holder Rod
- (1) TMM-1235-17 – Paper Roll Holder Mounting Bracket
- (4) TMM-1536-04 – 1/4-20 Hex Acorn Lock Nut
- (4) TMM-422-04 – Screw, 1/4-20 x .625 B.H.C.
- (2) TMM-681-04 – 1/4" Flat Washer

PREPARING THE BACKREST FOR THE PAPER ROLL HOLDER

Remove the back section pad from your TMM4 chair.

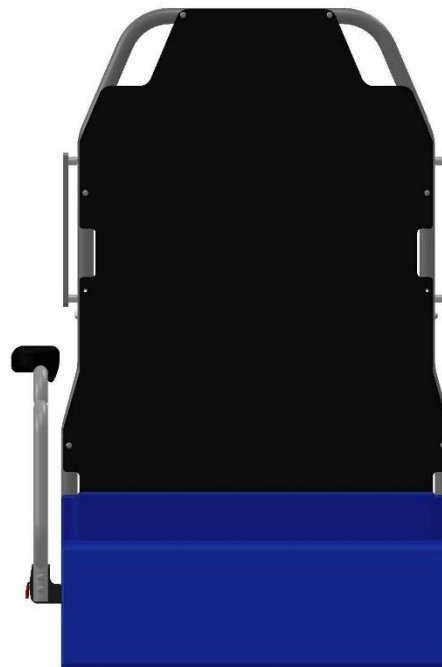

Raise the back section to the full upright position.

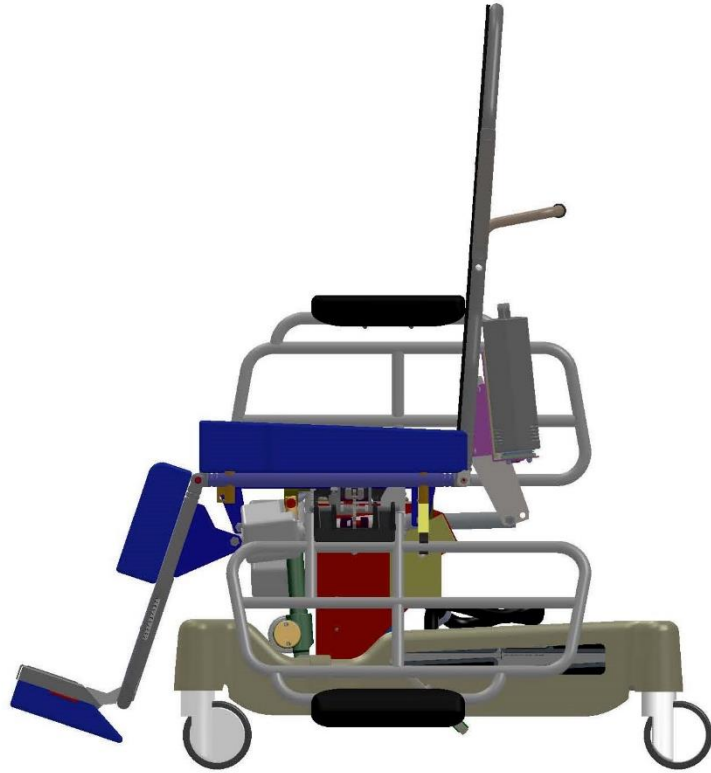

From the back of the chair, hold the Mounting Bracket (TMM-1235-17) against the phenolic board with the top edge of the bracket touching the back weldment tubing. Make sure the bracket is centered between the frame tubes and is parallel with the top of the chair.



Using the pencil, and the Paper Roll Mounting Bracket as a template, mark the two holes on the phenolic board and remove the bracket.

Using the 5/16" drill bit and Power Drill, drill two holes through the phenolic board.




NOTICE

Be careful to not force the drill through too quickly as you can chip the phenolic board.



INSTALLING THE PAPER ROLL HOLDER ASSEMBLY

Remove the Paper Roll Holder Rod (TMM-559-11) from the Paper Roll Bracket Assembly.

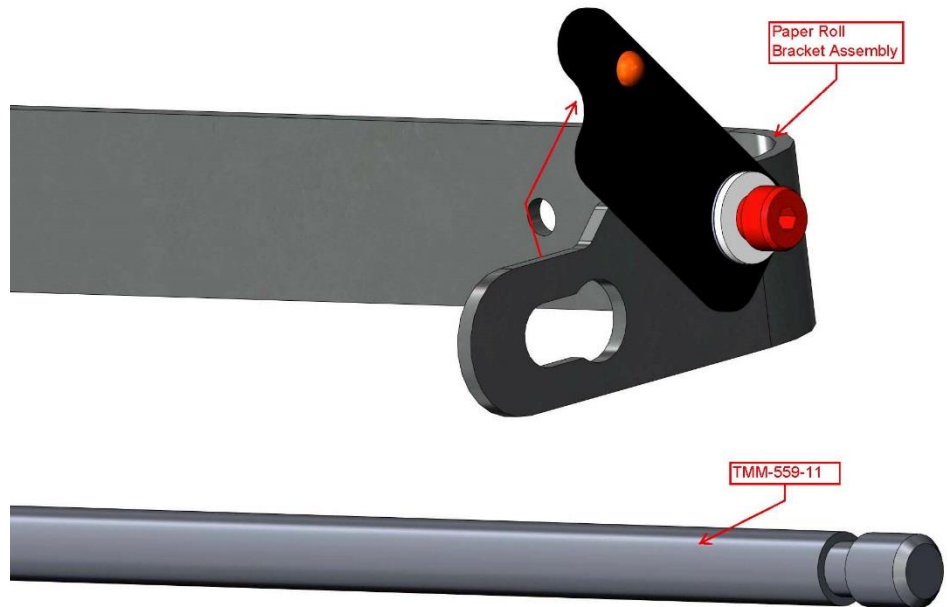


Using (2) TMM-422-04 screws, (2) TMM-681-04 flat washers and (2) TMM-1356-04 Acorn Nuts, mount the Paper Roll Bracket Assembly to the Paper Roll Holder Mounting Bracket (TMM-1235-17).



NOTICE

The Button Head Screws should be installed with the head on the inside of the Paper Roll Holder Bracket as shown.

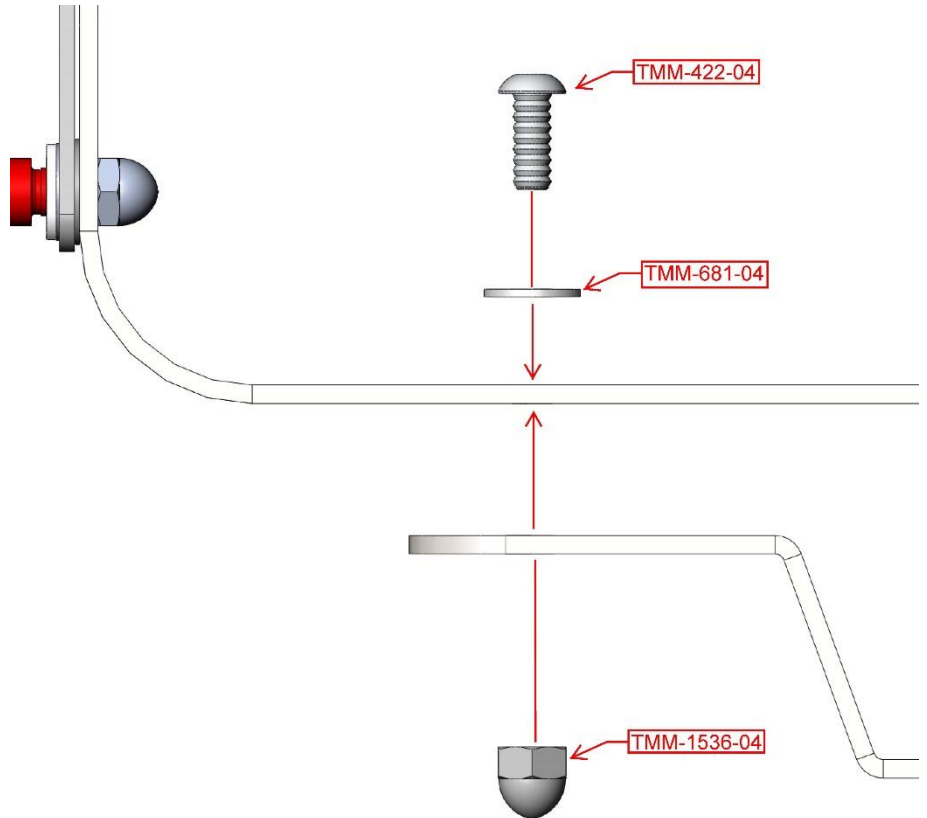

Using the Allen Wrench and 7/16" wrench, tighten the bolt and nuts securely.

Using (2) TMM-422-04 screws and (2) TMM-1536-04 Acorn Lock Nuts, install the Paper Roll Holder Assembly (TMA29-15) to the holes drilled in the phenolic board.




NOTICE

Be sure you have the Lock Levers on the Paper Roll Holder facing up towards the head end as shown.



NOTICE

The Button Head Screws should be installed with the head on the pad side of the board as shown.

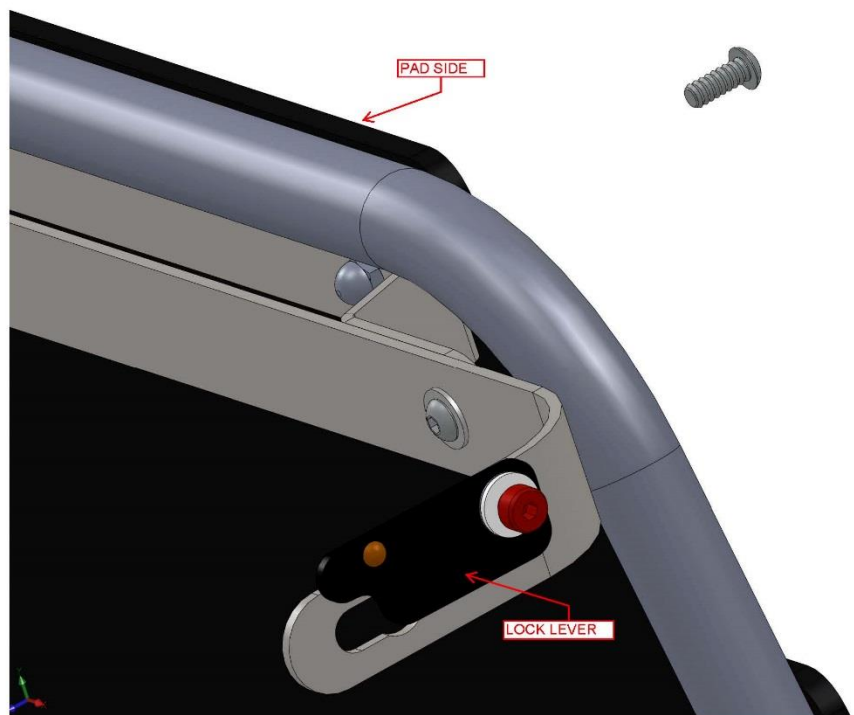

Using the Allen Wrench and 7/16" wrench, tighten the bolt and nuts securely.

CAUTION

Tighten fasteners periodically as routine maintenance of the accessory. Loose bolts may cause safety problems.

USING THE PAPER ROLL HOLDER ASSEMBLY


Lift the Lock Levers on both sides of the Paper Roll Holder.



Slide the Paper Roll Holder Rod into the large opening in the key slot and remove the rod.




Slide a roll of paper onto the rod.

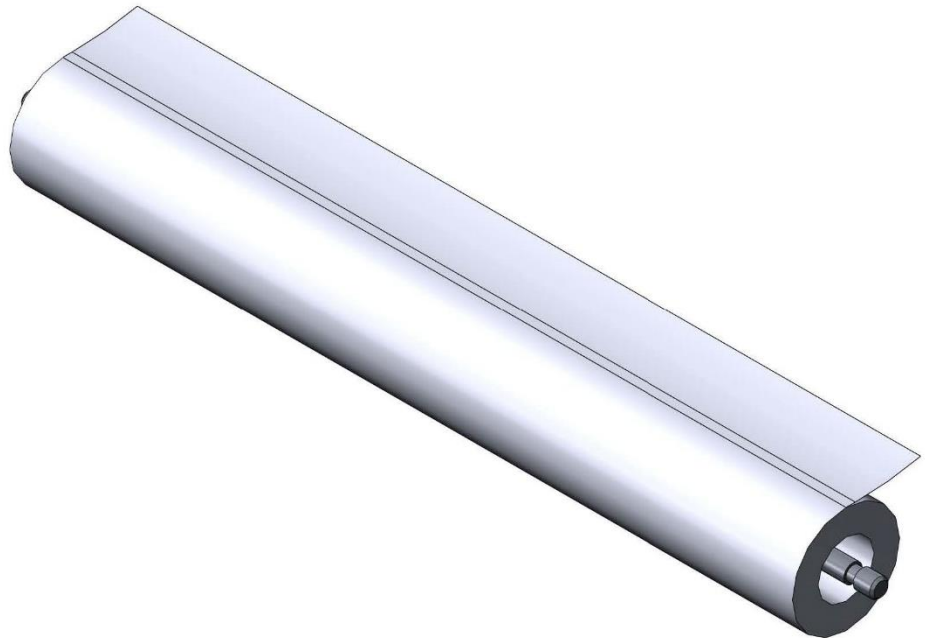

Feed one end of the rod into the large end of the keyhole shaped slot, sliding it past the groove and then insert the other end.




Center the rod and pull the rod back with the grooves into the smaller end of the keyhole shaped slot.

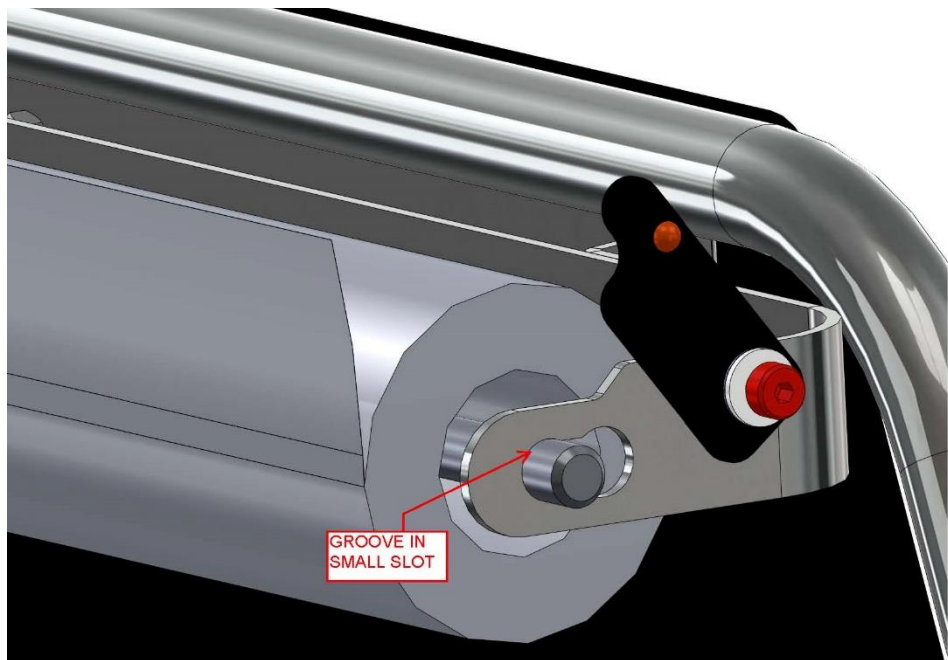

Lower both lock levers to prevent the rod from falling out.

CAUTION

If this accessory is exposed to bodily substances, clean thoroughly to prevent the spread of diseases.

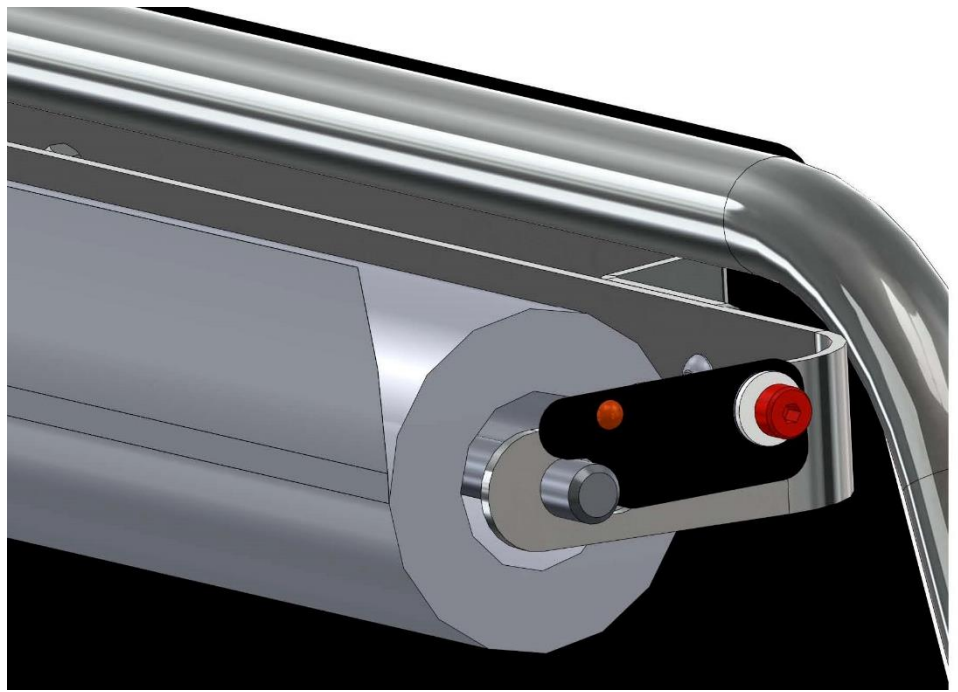

CLEANING INSTRUCTIONS



STEAM AND/OR PRESSURE CLEANING CHAIR WILL VOID WARRANTY

COMPONENT	FREQUENCY	PROCEDURE	CLEANING AGENT *	COMMENT
All Rails and Painted Metal	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Check for chipped paint/chrome
Base Cover	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Check for peeling or missing labels
Actuators and Battery Pack	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Clean exterior surfaces only with minimal water
Pendants	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Clean exterior surfaces only with minimal water
All Other Surfaces	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Check for damage
TransMotion Medical Accessories	As required/in accordance with facility policies	Wipe with damp cloth and allow to air dry	Facility approved detergents, disinfectants & water	Check for damage
3rd Party Accessories	As required/in accordance with facility policies	Consult manual provided by accessory manufacturer.		

* Use of Facility approved detergents and disinfectants must comply with the instructions provided by the manufacturer(s) of those products.

Upholstery Care & Cleaning Instructions

Proper care is essential in ensuring the durability and reliability of TransMotion Medical upholstery. In general, all products should be:

1. Cleaned
2. Disinfected (in accordance with facility policy)
3. Rinsed
4. Allowed to air dry

It is important to note that the terms “disinfecting” and “cleaning” should not be used interchangeably. Disinfectants alone will not provide adequate cleaning since they do not have the appropriate properties to cut grease or oil and remove grime, hair or skin oils. The appearance, feel and performance of your upholstery may diminish if not cleaned properly.

Cleaning:

Remove hair, grime and body oils. Your upholstery should be cleaned on a regular basis with a damp cloth soaked in a mild soap and water solution. Avoid harsh detergents or chemicals that could damage the finish of your product. If disinfecting with manufacturer approved chemicals or bleach, the upholstery MUST be wiped off using clean water on a damp cloth and allowed to air-dry as a final step. Failure to rinse upholstery with clean water can result in a build-up of residues that, over time, may lead to drying, cracking or other undesirable changes to appearance, feel and performance.

Retain all instructional tags for future use.

Contact the manufacturer of the fabric used in your upholstery for a current list of approved disinfectants. TransMotion Medical Customer Care can assist with identifying which upholstery was used with your Stretcher-Chair.

For customer supplied and non-standard materials, please refer to the individual manufacturer’s cleaning instructions.