

5560

Comfort Care, Premier Recliner



CUSTOMER INSTRUCTIONS

PLEASE READ AND FAMILIARIZE YOURSELF WITH ALL INSTRUCTIONS BEFORE USING THIS PRODUCT. If you have trouble understanding these instructions contact your dealer or Winco customer support, (800) 237-3377 before attempting to use this product; otherwise injury may occur.

Winco assumes no responsibility for damage or injury caused by the improper assembly, installation or use of these products; or during assembly or maintenance of these products.



1. READ AND FOLLOW ALL DIRECTIONS.
2. NEVER use side table as a seat.
3. DO NOT put hands, feet or clothing into any openings when changing positions on recliner. Attendant MUST confirm that users hands and feet are safe while changing recliner positions or INJURY MAY RESULT.
4. Stay clear of the recline mechanisms.
5. Periodically, recheck tightness of all nuts, bolts and screws.
6. NEVER use the trays or chair arms or backrest or legrest as a seat; SERIOUS INJURY MAY OCCUR.
7. Chair MUST be in full upright position with casters locked when a patient enters or exits chair.
8. NEVER allow a patient to exit a reclined chair with position-lock engaged or SERIOUS INJURY MAY RESULT.
9. Lock casters at all times, except when transporting chair.
10. Chair must ALWAYS be positioned on a level surface.
11. DO NOT use recliner for Transporting in or with ANY type of vehicle or trailer. Winco recliners have not been tested or approved for use by an occupant in any type of vehicle or trailer.
12. Immediately remove from service; Any recliner with broken recline mechanisms, torn upholstery, or other mechanical or visible damage.
13. USE ONLY WINCO AUTHORIZED REPLACEMENT PARTS.
14. NEVER EXCEED the recommended **weight capacity of 275 lbs.** (124.7 kg).
15. SAVE THESE INSTRUCTIONS for future reference and training.

5560 Weight Capacity = 275 lbs. (124.7 kg)

NO TOOLS REQUIRED FOR ASSEMBLY!

PLEASE READ AND FAMILIARIZE YOURSELF WITH ALL INSTRUCTIONS BEFORE PROCEEDING WITH ASSEMBLY

Install the Back:

1. After checking your product for any shipping damage, chair needs to be positioned upright on level surface. Remove all packing material and hardware that was secured for shipping. Remove tape and plastic. If using scissors or knife, please be careful not to damage upholstery.
2. Make sure footrest is in closed position & LOCK ALL CASTERS before installing chair back.
3. Remove the (2) **Bow Tie Cotter Pins** and the (2) **Clevis Pins** from the lower **Back Frame Forks**. **DO NOT LOOSE THESE PARTS, YOU WILL NEED THEM AGAIN!**
4. Using a second person to help; Position the **Back Frame Forks** over the **Seat Frame** and align holes.

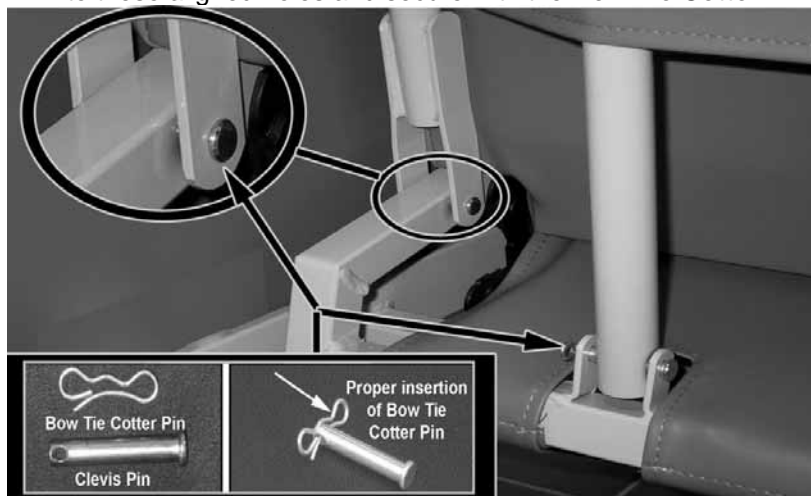
HINT: If you do not have a helper; Rest the bottom of the Upholstered Back on the rear portion of the Upholstered Seat and push the Back Assembly slightly forward to help steady the chair back.

5. Use (1) one of the **Clevis Pins** (with tapered end) and insert it THROUGH the inside of Back Frame Fork, THROUGH hole in Seat Rail Tube and; THROUGH outside of Back Frame Fork. REPEAT FOR OPPOSITE SIDE OF CHAIR. (see pg.2 detail of pin install photo)

**WARNING**

Failure to correctly install clevis pins connecting Back Frame fork holes and Seat Frame Rail holes, will cause back to operate poorly and may cause injury to the user or damage the chair and void warranty.

6. Insert (2) two **Bow Tie Cotter Pins**, into the small holes on ends of the (2) two **Clevis Pins** you just installed.
7. Remove the **Bow Tie Cotter Pin** and the **Clevis Pin** from the **Back Support Strut** (attached to backframe).
8. Align the holes of the **Back Support Strut** with the holes in the **Angle Brackets** (on chair frame) and insert the Clevis Pin into these aligned holes and secure with the **Bow Tie Cotter Pin**.

**Operation of the Chair:**

9. NON-RESTRAINT: To make the chair permanently non-restraint, remove the **Snap Button** of the **Right Tray Arm**, and unbolt the **Position Lock Bar** from the back of the chair.
10. This recliner has three (3) positions; UPRIGHT, SEMI-RECLINE, and RECLINE, which can be achieved by user or attendant. The position lock bar can only be engaged or disengaged by attendant. Not all notches on Position Lock Bar are used for this model chair.
11. POSITION LOCK: Attendant may hold the chair in a reclined position by engaging the Position Lock. To safely accomplish this, recline the chair to the 2nd or 3rd recline position and lower the position lock with your hand until the notches on the Position Lock bars engage with the Position Lock Pins.
To release Position Lock, raise the Position Lock off of the Pins and return chair to upright position.

WARNING: Attendant must release Position Lock before patient can exit chair – If patient attempts to exit chair with Position Lock engaged, **SERIOUS INJURY MAY RESULT.**

WARNING: Failure to release the Position Lock before changing positions will result in damage to the chair and void the warranty.

12. This chair is equipped with **Total Lock casters**. Stepping on the brake tab will lock the wheel and the swivel mechanism of the caster for added safety and security.

To ensure proper operation of the reclining mechanism, lubricate all pivot points with a silicone type lubricant (or lubricant approved by your facility) at least twice yearly. Periodically recheck tightness of all screws and bolts.

Wrinkles in Vinyl – Storage and shipping can cause wrinkles of the vinyl on your Winco chair. This is normal. Please allow 1-2 weeks at room temperature (approx. 70 degrees) to allow for vinyl to expand back to original appearance

Periodically recheck tightness of all screws and bolts, and never use the chair arms, backrest, legrest or tray as a seat.

Visit our website www.wincomfg.com/product-documents/parts-lists for a complete parts list catalog of your chair or call Customer Service @ 1-800-237-3377 to request it by mail.



LIMITED WARRANTY

5516 SW 1st Lane • Ocala, FL 34474
Phone (352) 854-2929 • (800) 237-3377 • Fax (352)
854-9544

TERMS

Winco Mfg., LLC. warrants this product to be free of manufacturer's defects in material and workmanship, provided that the product is used according to normal operating conditions and proper maintenance intended by manufacturer. This warranty is available only to the original retail purchaser, is non-transferable and commences on the date of retail sale; proof of purchase required.

Limited One Year Complete Product Warranty: Winco Mfg., LLC warrants the complete product for one (1) year. At Winco Mfg., LLC. sole discretion, it may repair or replace any components freight free that are found to be defective during the first year. *Winco Mfg., LLC. shall not be liable for any labor, or any other costs incurred as a result of or in conjunction with a warranty claim.*

Limited Three Year Warranty: Winco Mfg., LLC warrants the recliner mechanisms, electrical components, vinyl [from cracking or delaminating] pressurized gas springs, when new, for a period of three (3) years. At Winco Mfg., LLC. sole discretion, it may repair or replace components found to be defective. *Winco Mfg., LLC. shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.*

Limited Lifetime Chair Frame Warranty: Winco Mfg., LLC warrants the steel chair frame, when new, for the lifetime of the original purchaser. This limited warranty does not apply to paint/finish or any components attached to the frame such as; upholstery, foam, casters, mechanisms or related parts that are covered under above warranties. *Winco Mfg., LLC shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.*

The purchaser's exclusive remedy under this warranty shall be limited to such repair or replacement of defective components at Winco Mfg., LLC sole discretion. For warranty service, contact the dealer that the product was originally purchased from or Winco Mfg., LLC directly.

EXCLUSIONS

There are no other warranties, conditions, representations or guarantees, express or implied, made or intended by Winco Mfg., LLC and all other warranties, conditions, representations or guarantees including any warranties, conditions, representations or guarantees under any Sale of Good Act or Like legislation or statute is hereby expressly excluded. Any and all other implied warranties shall not extend beyond the duration of the express warranty. Liability for incidental or consequential damages is excluded to the extent permitted by law. Some states do not allow the exclusion or limitation so the above limitation or exclusion may not apply to you. This warranty gives you specific legal rights. In addition you may also have other rights, which vary from state to state.

GENERAL PROVISIONS

No warranty herein contained or set out shall apply when damage or repair is caused by any of the following:

1. Damage in transit.
2. Accident, alteration, abuse or misuse of product ie. exceeding weight capacities, applying inordinate pressure to footrest/leg rest, use in Mental Health facilities or like institutions.
3. Fire, water damage, theft, war, riot, hostility, acts of God.

Examples contained in this list are not to be construed as all-inclusive.

